# EVALUATION OF FECAL PYRUVATEKINASE ISOENZYME (M2-PK) LEVEL IN DIFFERENTIATING FUNCTIONAL FROM ORGANIC COLONIC DISORDERS

By

# Abd Al-Azim Hassan<sup>1</sup>, Mohammed Bastawy<sup>2</sup>, Moustafa El-Hawary<sup>3</sup> and Kamel Soliman<sup>4</sup>

Departments of Tropical Medicine and Clinical Pathology, Faculty of Medicine, Al-Azhar University

# ABSTRACT

**Background:** Pyruvatekinase, which converts phosphoenolpyruvate to pyruvate, is a key enzyme in glucose metabolism and is present in organ-specific isoforms (the L, R, M1, and M2 isoforms). In normal proliferating cells, M2-PK is mainly tetrameric and has a high affinity for phosphoenolpyruvate. In contrast, the M2-PK isoenzyme found in tumor cells is usually dimeric and has a low affinity for phosphoenolpyruvate. For this reason, the dimeric form of M2-PK has been named tumor M2-PK.

**Objective:** This study evaluated the potential value of fecal, dimeric M2-PK level in differentiating functional from organic colonic disorders as well as its value as a surrogate marker of inflammation in patients with inflammatory bowel disease (IBD) and colorectal cancer (CRC).

**Patients and Methods:** This prospective study included 60 patients with different colonic disorders, 20 patients with Functional colonic disorders, 20 patients with inflammatory bowel disease (ulcerative colitis & chron's disease), and 20 patients with colon cancer. The M2-PK level was measured in all patients using a highly sensitive enzyme – linked immunosorbent assay (ELISA), which allowed the quantitative measurement of tumor M2-PK in stool.

**Results:** Our study revealed a highly significant increase in tumor M2-PK in the stool samples of those patients with organic colonic disorders (IBD and CRC groups) compared to functional group (IBS). At a cutoff value of 4.2 (U/ml), our overall sensitivity and specificity for organic group over the functional group were 87.5% and 80% respectively. Furthermore, the results of M2-PK levels (U/mL) in our study were shown to be significantly elevated in active, compared to inactive IBD.

**Conclusion:** In colonic disorders, fecal concentration of tumor M2-PK was a good marker for discrimination of functional from organic colonic conditions (IBD and CRC), with a sensitivity and specificity of **87.5% and 80%** respectively. Tumor M2-PK can be used as a tumor marker in screening of colorectal cancer. Dimeric fecal M2-PK has the potential to be an important noninvasive marker of disease activity in IBD.

Key words: M2-pyruvate kinase, inflammatory bowel disease, Crohn's disease, ulcerative colitis, colorectal cancer, irritable bowel syndrome.

## **INTRODUCTION**

Differentiating between inflammatory bowel disease (IBD) and functional bowel

disorders, such as irritable bowel syndrome (IBS), can often be difficult as they present with similar symptoms *(Kenneth, 2011).* 

The diagnosis of IBD typically necessitates invasive endoscopic procedures to visualize the mucosa and enable confirmatory histological specimens to be obtained. However, this may miss disease in the gastrointestinal (GI) tract, not directly visualized endoscopy. at Furthermore, in IBD accurate monitoring of disease activity may include repeated endoscopy, as symptoms correlate poorly with disease activity (Sutherland et al., 2008).

Noninvasive biomarkers in IBD are being increasingly recognized as important, both at the initial diagnosis and for monitoring disease activity. They also play a valuable role in differentiating organic GI disease from functional disorders by examining the entire GI tract (*D*?britz et al., 2014).

Key characteristics of fecal biomarkers include stability in fecal samples and the existence of a sensitive and reliable assay (*Christofk et al., 2008*).

Pyruvatekinase, which converts phosphoenolpyruvate to pyruvate, is a key enzyme in glucose metabolism and is present in organ-specific isoforms (the L, R, M1, and M2 isoforms). In normal proliferating cells, M2-PK is mainly tetrameric and has a high affinity for phosphoenolpyruvate. In contrast, the M2PK isoenzyme found in tumor cells is usually dimeric and has a low affinity for phosphoenolpyruvate. Dissociation of the tetrameric form to the dimeric form in tumor cells is induced by direct interaction of M2-PK with various oncoproteins. For this reason, the dimeric form of M2PK has been named tumor M2PK. Because of its low affinity for phosphoenolpyruvate, tumor M2-PK is easily released from tumor cells and is quantitatively detectable in body fluids (*Abdullah et al.*, 2012).

Tumor M2-PK can also be detected and quantified in stool samples using an ELISA. It was shown that fecal tumor M2PK is more accurate for CRC screening than serum or plasma tumor M2PK (*Tonus et al., 2012*). The fecal tumor M2-PK test has good sensitivity and specificity for CRC detection, especially in high-risk or symptomatic populations (*Abdullah et al., 2012*).

M2-PK is crucial for rapid tumor growth and aerobic glycolysis during tumorigenesis (*Christofk et al., 2008*). Its tetrameric and dimeric forms are present in rapidly proliferating cells of many tissues, including leukocytes. Upon leukocyte destruction in the gastrointestinal tract, the protein is released to the fecal stream (*Gupta and Bamezai*, 2010).

M2-PK is stable in stools, which increases the potential value of fecal M2-PK concentration assessment both in intestinal inflammation and cancer. They postulated that fecal M2-PK could serve as a biomarker of inflammation in inflammatory bowel diseases (IBD) (*Czub et al., 2007*).

Active IBD is accompanied by increased cell turnover and rapid division (*Foell et al., 2009*). *Turner et al. (2010*) investigated fecal M2-PK, FC, lactoferrin and S100A12 protein in children with severe ulcerative colitis. They showed that only M2-PK had constructive and predictive validity, while other markers failed to meet this criterion. The present work aimed to evaluate the potential value of fecal, dimeric M2-PK level in differentiating functional from organic colonic disorders, as well as its value as a surrogate marker of inflammation in patients with inflammatory bowel disease (IBD) and colorectal cancer (CRC).

# **PATIENTS AND METHODS**

Patients attending the gastroenterology and outpatient clinic in patient Department of Gastroenterology and Oncosurgery, El Saved Galal teaching hospital, during the period from August 2014 to August 2016, with new-onset lower GI symptoms, with previously diagnosed with previously IBD or prospectively diagnosed CRC were enrolled in the study. Patients with dyspepsia symptoms of or gastroesophageal reflux disease were excluded. All patients underwent clinical evaluation, including history (including medication) and full physical examination before providing a fecal sample for M2-PK assays. Patients presenting with diarrhea also provided stool samples for microscopy (ova, cysts, and parasites). Further investigations, including endoscopic examinations, were enrolled as clinically indicated. For the purposes of this analysis, patients were subdivided into the following diagnostic categories: IBD group (ulcerative colitis (UC) and Crohn's disease (CD), irritable bowel syndrome/functional bowel disorder (IBS), and colorectal cancer (CA). In a subanalysis of IBD patients, disease activity was assessed by a combination of physician global assessment, and/or endoscopic grading where available (D'Haens et al., 2007). Histopathological

examinations of colonic biopsies were done.

Patients provided a single stool sample for analysis, submitted within 48 hours. For M2-PK, a commercially available ELISA kit (ScheboTech, Giessen. Germany), which has no cross reactivity to other forms of pyruvate kinase was used. Patients did not have to keep a special diet and were told to take their usual medications. Usually 10mg feces samples were collected by sterile tube. This assay has an intratest and interest variability coefficient of 4.5% and 6.1%, respectively. Stool samples for M2-PK are stable at room temperature (21°C) for 3 days or for up to 1 year at -20°C (Tonus et al., 2012).

This study was approved by Al-Azhar Faculty of Medicine Local Research Ethics Committee. The protocol was explained to patients selected for the study or at least one of their relative, and written informed consents were obtained from all patients who participated or from their relatives.

**Statistical analysis:** The data were analyzed using Microsoft Excel 2010 and statistical package for social science (SPSS version 24.0) for windows (SPSS IBM., Chicago, IL).

Results were expressed as mean  $\pm$  SE with 95% confidence interval using mean for quantitative variables, and using the frequencies and percentage for qualitative ones; a P value < 0.05 was considered statistically significant.

Quantitative data were analyzed by applying one way analysis of variance (ANOVA) test for comparison of the mean of more than two groups, while Student t-test was used for comparison of the mean of two groups.

The ROC curve was used to detect the cutoff points and the correlation between the sensitivity and specificity of functional, IBD and CRC groups.

### RESULTS

Of the 60 patients evaluated, 18 had biopsy-proven UC, 2 had biopsy-proven CD and 20 had biopsy-proven CRC. In all, 20 had normal investigations and were diagnosed as having functional bowel disorders/IBS after full assessment.

In the functional bowel disorder subgroup, the vast majority of patients had IBS according to ROME III criteria. However, in all these patients all clinically indicated endoscopic investigations were normal. Student t-test was positive for age between IBS group and each one of the other diagnostic groups, indicating a significant difference.

Table (1): Age and gender in the Main Diagnostic Groups

Groups	Ν	M:F Ratio	Age (yrs)	P for age* (vs. IBS)
Colorectal cancer	20	14:6	40-79	0.03
Inflammatory bowel	20	6:14	18 – 57	0.03
disease				
<b>IBS/functional bowel</b>	20	10:10	21 - 55	-
Total	60	30:30		-

Table (2): Fecal M2PK Values in the Main Diagnostic Groups:

Fecal M2PK						P. Value		
(U/mL)	Ν	Min	Max	Mean	S.E	a	b	
Groups								
<b>IBS/functional bowel</b>	20	1.20	29.0	5.18	1.5			
Inflammatory bowel disease	20	3.0	279.0	49.9	15.5	0.01*		
Colorectal cancer	20	1.6	217.0	62.5	14.3	0.001**	0.5	
<sup>a</sup> p= Another groups compared with Functional group.								
<sup>b</sup> p= CRC group compared with IBD group.								

Kruskal–Wallis ANOVA was positive for M2-PK between the 3 diagnostic groups, indicating a significant difference between the groups. The mean values of M2-PK were highly significantly elevated in IBD and CRC patients compared to IBS. Although M2-PK concentrations in IBD patients were higher than the CRC group; these differences were not statistically significant (**Table 2**).

There were highly significant differences in the M2-PK concentrations in patients with active versus inactive diseases (Table 3).

Fecal M2PK (U/mL) IBD	N	Min	Max	Mean	S.E	P. Value
In remission	10	7.0	35.0	16.7	2.6	0.02*
In exacerbation	10	3.0	279.0	83.3	27.6	0.02

Table (3): Fecal M2PK Values in Active vs. Inactive IBD disease:

Fecal M2-PK had a sensitivity of 87.5%, specificity of 80%, PPV of 89.7%, and a negative predictive value (NPV) of 76.2% for organic GI diseases. Similar results were obtained when the analyses were performed for IBD alone (sensitivity of 95%, specificity of 80%, PPV of 82%,

and a negative predictive value (NPV) of 94%). For colorectal cancer patients alone, sensitivity of 87.5%, specificity of 80%, PPV of 89.7%, and a negative predictive value (NPV) of 76.2% (**Fig. 1** and Table 4).

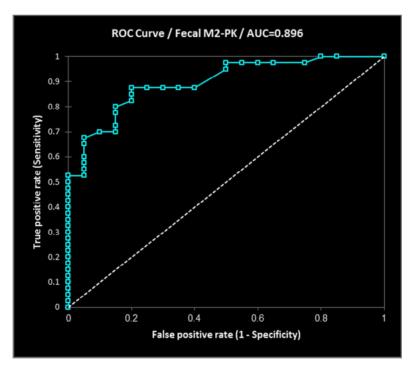


Figure (1): Receiver operator characteristic (ROC) plot for Fecal M2-PK between Functional and Organic groups (area under the curve AUC = 0.896)

**Table (4):** Sensitivity, specificity, efficacy, positive and negative predictive values of the plasma concentrations of tumer M2-PK for Organic GI Disease.

Cutoff	Sensitivity	Specificity	PPV	NPV	ТР	TN	FP	FN	Accuracy	P. Value
4.2	87.5	80.0	89.7	76.2	35	16	4	5	85.0	0.001**

#### **DISCUSSION**

The gold standard investigation for the early detection of colorectal cancer is colonoscopy. However, the acceptance of this costly and invasive method is low (*Mariann et al., 2016*).

In order to increase the participation in colorectal cancer screening programs, an easy, fast and economical initial screening method, with good patient compliance, is absolutely necessary. This allows identification of those patients most likely to have colorectal cancer, who require further investigation by colonoscopy.

An increase in tumor M2-PK in EDTA plasma samples is found in gastrointestinal cancers, as well as a wide range of other tumors such as lung, renal, breast and cervical cancer. Thus, fecal tumor M2PK is more accurate for CRC screening than serum or plasma tumor M2 PK (*Tonus et al., 2012*).

Our study revealed a highly significant increase in tumor M2-PK in the stool samples of those patients with organic colonic disorders (IBD and CRC groups) compared to functional group (IBS). At a cut-off value of 4.2 (U/ml), our overall sensitivity and specificity for organic group over the functional group were 87.5% and 80% respectively. These data corresponded well with Chung-Faye et al. (2007) who found that mean M2-PK values were significantly elevated in UC, CD, and colorectal carcinoma compared to IBS. Using a predetermined cut-off level for normal fecal M2-PK. a sensitivity, specificity, and positive predictive value of 73%, 74%, and 89%, respectively for differentiating organic disease from IBS were obtained.

These data also corresponded well with *Mariann et al.*, (2016) who found that sensitivity and specificity of M2-PK for CRC were 94.7% and 67.5% respectively.

Studies suggest that inflammatory reactions in the bowel can cause an elevation in fecal tumor M2 PK level (*Mulder et al., 2007 and Shastri et al., 2008*). In our study, 95% of patients with IBD had positive M2 PK test results at cutoff value of **4.2** U/ml. Most IBD patients with positive results had ulcerative colitis IBD.

The results of the present study showed significant differences in the level of fecal tumor M2-PK between IBD cases and functional group (IBS) cases, while the results revealed non-significant а difference between IBD and CRC groups. Of the 20 patients evaluated with IBD, 18 had biopsy-proven ulcerative colitis and 2 had biopsy-proven crohn's disease. The mean tumor M2-PK value in IBD group was 49.9 U/ul, and at a cut-off value of 4.2 U/ml. Our overall sensitivity for IBD was 95 %, and its specificity was 80 %. These data came in correlation with the results of a study done by Mulder et al. (2007) who showed that in IBD patients, the M2-PK test was positive in (78.9%). In another multicenter study done by Shastri et al. (2008), up to 88% of patients with ulcerative colitis and Crohn's Disease had elevated M2-PK levels. Furthermore, the results of M2-PK levels (U/mL) in our study were shown to significantly elevated in active, be compared to inactive, IBD and these data came in correlation with Chung-Faye et al. (2007).

In our study, the functional group consisted of 20 individuals with ROME

III citeria for IBS. The mean tumor M2-PK value in this group was 5.18 U/ml. In 16 of 20 subjects, tumor M2-PK levels were below the cut-off value. The resulting specificity at a cut off value of 4.2 U/ml was 80 %, which is in general accordance with the study of Koss et al. (2008) who reported specificities between 50% and 98%. Positive results in patients with normal colonoscopy results may have been caused by inflammation, which have been detected by microscopic examination. Shastri et al. (2008) also showed that 26.2% of normal subjects, including patients with irritable bowel syndrome, (135/156; cutoff value >4 U/mL) had positive fecal tumor M2-PK test results. By using a combination of clinical factors, such as age, symptoms (especially those fulfilling ROME II criteria), and noninvasive fecal markers, a significant proportion of patients with inactive IBD or IBS may avoid further endoscopic investigations, invasive reserving endoscopy only for those with elevated fecal M2-PK concentrations.

Cancer cells have higher rates of glucose cells. uptake than normal However, only a small fraction of the glucose taken up is used for oxidative phosphorylation. The decrease in aerobic glycolysis in these cells may be due to reprogramming of metabolic genes, enabling cancer cells to partition glucose metabolites away from oxidation towards the synthesis of macromolecules (Uppara et al., 2015). M2-PK is crucial for rapid tumor growth and aerobic glycolysis during tumorigenesis (Christofk et al., 2008). Therefore, M2-PK can be used as a noninvasive biomarker and a diagnostic tool for the detection of cancer.

## CONCLUSION

Fecal concentration of tumor M2-PK was a good marker for discrimination of functional from organic colonic conditions (IBD and CRC); with a sensitivity and specificity of **87.5% and 80%** respectively. Tumor M2-PK can be used as a tumor marker in screening of colorectal cancer. Dimeric fecal M2-PK has the potential to be an important noninvasive marker of disease activity in IBD.

#### REFERENCES

- 1. Abdullah M, Rani AA, Simadibrata M, Fauzi A and Syam AF (2012): The value of fecal tumor M2 pyruvate kinase as a diagnostic tool for colorectal cancer screening. Acta Med Indones. 44(2):94-9.
- Christofk HR, Vander Heiden MG, Harris MH, Ramanathan A, Gerszten RE, Wei R, Fleming MD, Schreiber SL and Cantley LC (2008): The M2 splice isoform of pyruvate kinase is important for cancer metabolism and tumour growth. Nature, 452: 230-233.
- **3.** Chung-Faye G, Hayee B, Maestranzi S, Donaldson N, Forgacs I and Sherwood R (2007): Fecal M2-pyruvate kinase (M2-PK): a novel marker of intestinal inflammation. Inflamm Bowel Dis., 13:1374–1378.
- 4. Czub E, Herzig KH, Szaflarska-Popawska A, Kiehne K, Socha P, Woś H, Kamińska B, Błaszczyński M, Cichy W, Bała G, Brodzicki J, Grzybowska U and Walkowiak J (2007): Fecal pyruvate kinase: a potential new marker for intestinal inflammation in children with inflammatory bowel disease. Scand J Gastroenterol., 42: 1147-1150.
- **5. D?britz J, Musci J and Foell D (2014):** Diagnostic utility of faecal biomarkers in patients with irritable bowel syndrome World J Gastroenterol., 20(2): 363-375.
- D'Haens G, Sandborn WJ, Feagan BG, Geboes K, Hanauer SB, Irvine EJ, Lémann M, Marteau P, Rutgeerts P, Sch?lmerich J and Sutherland LR (2007): A review of

activity indices and efficacy end points for clinical trials of medical therapy in adults with ulcerative colitis. Gastroenterology, 132:763-86.

- **7. Foell D, Wittkowski H and Roth J (2009):** Monitoring disease activity by stool analyses: from occult blood to molecular markers of intestinal inflammation and damage. Gut, 58: 859-868.
- **8.** Gupta V and Bamezai RNK (2010): Human pyruvate kinase M2: a multifunctional protein. Protein Sci Publ Protein Soc., **19**: 2031–2044.
- **9. Kenneth R (2011):** Gastrointestinal Disorders. In: Current diagnosis and treatment. 50th ed.pbi. Lange medical books: McGraw Hill, USA; PP. 609-10.
- **10.** Koss K, Maxton D and Jankowski JA (2008): Faecal dimeric M2 pyruvate kinase in colorectal cancer and polyps correlates with tumour staging and surgical intervention. Colorectal Dis., 10:244–248.
- 11. Mariann R, Ren?ta B, Anita B, Anna F, ?gnes M, Ferenc N, Zolt?n S, M?nika S, L?szl? T, Klaudia F and Tam?s M (2016): Diagnostic Accuracy of Five Different Fecal Markers for the Detection of Precancerous and Cancerous Lesions of the Colorectum. Mediators of Inflammation. Article ID 2492081, 6 pages.
- 12. Mulder S, Van ME, Anneke J, Jan F, Albert W, Ernst J and Rob J (2007): Tumor pyruvate kinase isoenzyme type M2 and immunochemical fecal occult blood test: performance in screening for colorectal cancer. Eur J Gastroenterol Hepatol., 19:878–82.

- **13. Shastri YM, Povse N, Schroder O and Stein J (2008):** Comparison of a novel fecal marker—fecal tumor pyruvate kinase type M2 (M2-PK) with fecal calprotectin in patients with inflammatory bowel disease: a prospective study. Clin Lab., 54:389–90.
- **14. Sutherland AD, Gearry RB and Frizelle FA** (2008): Review of fecal biomarkers in inflammatory bowel disease. Dis. Colon Rectum., 51:1283–91.
- **15.** Tonus C, Sellinger M, Koss K and Neupert G (2012): Faecal pyruvate kinase isoenzyme type M2 for colorectal cancer screening: a meta-analysis. World J Gastrointest Oncol., 18(30):4004–11.
- 16. Turner D, Leach ST, Mack D, Uusoue K, McLernon R, Hyams J, Leleiko N, Walters TD, Crandall W, Markowitz J, Otley AR, Griffiths AM and Day AS (2010): Faecal calprotectin, lactoferrin, M2-pyruvate kinase and S100A12 in severe ulcerative colitis: a prospective multicentre comparison of predicting outcomes and monitoring response. Gut, 59: 1207-1212.
- 17. Uppara1 M, Adaba F, Askari1 A, Clark S, Hanna G, Athanasiou T and Faiz1 O (2015): A systematic review and meta-analysis of the diagnostic accuracy of pyruvate kinase M2 isoenzymatic assay in diagnosing colorectal cancer. World Journal of Surgical Oncology, 13:48.

عبد العظيم حسن - محمد بسطاوي- مصطفى الهوارى-كامل سليمان

قسمى الأمراض المتوطنة والباثولوجيا الإكلينيكية- كلية طب الأزهر

خلفية البحث : يعتبر البيروفات كيناز اللذى يقوم بتحويل الفوسفو إينول بيروفات إلى بيروفات إنزيم رئيسى فى عملية حرق الجلوكوز ويوجد فى أشكال مرتبطة بنوع العضو (ل، ر، م ۱ و م ۲). فى الخلايا ذات التكاثر الطبيعى يكون مشتق البيروفات كيناز (م ۲) على شكل ثلاثى وله انجذاب شديد للفوسفو اينول بيروفات. على النقيض يوجد مشتق البيروفات كيناز (م ۲) فى الخلايا الورمية على شكل ثنائى وله إنجذاب ضعيف للفوسفو إينول بيروفات، ولذلك تمت تسميته البيروفات كيناز م٢

**الهدف من البحث:** تقييم معدل مشتق البيروفات كيناز في البراز في التمييز بين أمراض القولون الوظيفية والعضوية ودوره كمؤشر بديل للإلتهاب في مرضى إلتهابات الأمعاء وأورام القولون والمستقيم

المرضى وطرق البحث : شملت هذه الدراسة المستقبلية ٢٠ مريضا يعانون من أمراض قولون متنوعة. ٢٠ مريضا يعانون من أمراض قولون وظييفية، ٢٠ مريضا يعانون من أمراض إلتهابات الأمعاء (القولون التقرحى ومرض كرونز)، و ٢٠ مريضا مصابين بسرطان القولون. وقد تم قياس معدل البيروفات كيناز (م٢) لكل المرضى باستخدام التقنية عالية الحساسية (الإليزا) والتي تسمح بالقياس الكمي للبيروفات كيناز (م٢) في البراز.

النتائج: أظهرت هذه الدراسة وجود زيادة عالية ذات دلالة لمشتق البيروفات كيناز فى عينات براز المرضى المصابين بأمراض القولون العضوية (أورام القولون والتهابات القولون المناعية)، مقارنة بمجموعة أمراض القولون الوظيفية (متلازمة الأمعاء العصبية)، وكانت الحساسية و النوعية لمشتق البيروفات كيناز لانتقاء المرضى المصابين بأمراض القولون العضوية ٥٩٨ % و ٨٠%، عند قيمة ٤,٢ وحدة / مللى لتر. كما كانت نتائج معدلات البيروفات كيناز فى دراستنا هذه عالية وذات دلالة احصائية فى مجموعة المرضى المصابين بإلتهابات القولون المناعى النشط مقارنة بالمرضى المصابين بإلتهابات القولون المناعى الخامل.

الخلاصة : مستويات البيروفات كيناز فى البراز يمكن أن تكون ذات دلالة لتمييز أمراض القولون الوظيفية من أمراض القولون العضوية (إلتهابات القولون المناعية وأورام القولون والمستقيم)، بمدى حساسية ونوعية ٨٧,٥% و ٨٠% على التوالى. ويمكن استخدام مشتق البيروفات كيناز (٢٥) الورمى كدلالة لإستكشاف أورام القولون والمستقيم وقد بينت هذه الدراسة أن مشتق البيروفات كيناز (٢٥) الثنائى له الإمكانية أن يكون علامة غير مؤذية لنشاط مرض إلتهابات الأمعاء.